

Receiver DH-10

Perfect for hearing impaired educationalists



# DH-10 Neck-Worn Receiver

The DH-10 Receiver can be used with all Comfort Audio Digisystem microphones for the clearest possible reception of speech.

The lightweight DH-10 can be used with either a neckloop and T-Coil enabled hearing aid, or with headphones or earbuds. Simple controls allow adjustment of volume and sound balance.

When used in conjunction with either the Comfort Audio Conference Microphone DC-10 or Microphone DM-20, the DH-10 provides the best listening solution for hearing impaired teachers and lecturers.

The microphone is an advanced microphone, sound processor and transmitter. It captures voices from your classroom, cleans it of all background racket and emphasises the speech frequencies, which it then transmits to your wireless receiver without any perceptible delay. Use the DC-10 in your classroom together with one of our digital receivers and you'll be amazed at how clearly you will hear what's being spoken.



## Unique characteristics

**SST [SecureStream Technology] – Digital sound transfer**  
This unique digital sound transfer between the transmitter and the receiver eliminates all transmission noise. Moreover, the signals are encrypted, so no one can overhear what is being said. Unlike other conventional encrypted systems, SST is free from time lag and crystal clear.

**FFC – Free Frequency Check**  
Indicates the suitable frequencies in the specific environment.

## Receiver DH-10

### General Technical Specifications

#### General

**Size:** 74 x 39 x 13 mm

**Weight:** 35 g

#### Special functions

**Balance:** Easily adjustable

**Fix volume:** Option to set fixed volume  
Volume control is switched off.

**Max volume:** Option to set max volume  
Volume control is limited.

**PIN-code:** Functions can be locked with  
a PIN-code.

#### User interface

**Display:** Coloured

**Menu:** Easy to navigate; automatic  
energy saver

**Menu languages:** English, Swedish,  
Danish, Norwegian, French, German,  
Dutch

**Function indicator:** Easy-to-read, coloured  
status information

**Manual:** USB or [www.comfortaudio.com](http://www.comfortaudio.com)

#### Power supply

**Battery:** Chargeable Li-Polymer

**Operating time:** Up to 10 hours

**Charging:** Overcharging protection,  
automatic float charging

**Charging time:** Approx. 2 hours, 90%  
after 1 hour; can be charged during use.

#### Radio Transmission

**SST:** Secure Stream Technology – FSK,  
encrypted digital sound transfer.

**FFC:** Free Frequency Check – indicates  
the suitable frequencies in the specific  
environment.

**Frequencies receiver:** 768 – 872 MHz

**Range:** Approx. 30 m

#### Technical specifications

**Type of receiver:** SST





**Antenna:** Built-in

**Connections:** Neck loop, headphones/  
earphones

#### Functions



#### Applications:

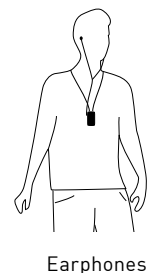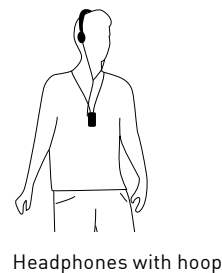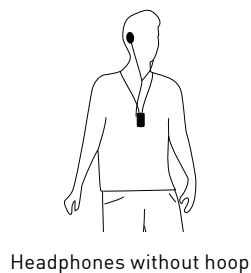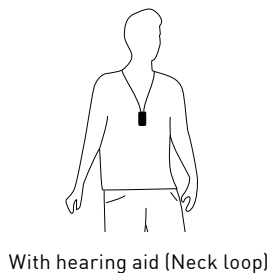

100224

SPECIFICATION CAN BE SUBJECT TO CHANGE WITHOUT NOTICE

THE UK'S LEADING SOURCE OF AUDIOLOGICAL INSTRUMENTS